

SWIPPER™

TESTING KITS
FOR DRUGS

Fentanyl Test Kits

SW-0690



**Read
Before
Use**

As a precaution, gloves and eye protection should be worn when handling substances and using test units. Do not puncture any part of unit. Do not remove probe after breaking the ampule. In event of a spill, rinse affected area thoroughly with water. Keep out of reach of children. Store the tubes in a dark place, not exposed to light. All results are presumptive. If necessary, further testing by a competent drug testing lab can provide additional identification.

Test Preparation

1

Remove the SWIPPER™ Cap from the test unit.




Ensure that the test tube remains upright in order to avoid the ampule falling out of the tube.

2

Sample the suspected substance with the SWIPPER™ Cap.


Sample multiple locations

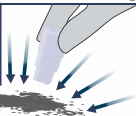

Test Preparation

3

Return the SWIPPER™ Cap to the test tube. Press firmly.

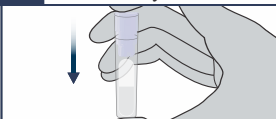

STOP

Do not continue unless the SWIPPER™ Cap is properly secured.

⚠ CAUTION ⚠

The kit contains an acid. Make sure to follow all the safety procedures.

Identification


4

With **thumb and forefinger**, squeeze the test tube until the ampule breaks.



5

Gently shake test up and down to ensure that the liquid touches the Cap.

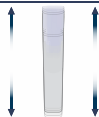

Results



The appearance of a red-purple (wine) color indicates the presence of **Fentanyl**.



The appearance of a reddish-orange color within one minute indicates the presence of **Amphetamine** or **Methamphetamine**.

* This does not exclude the possibility that Fentanyl may be present.



The appearance of a purple color within one minute indicates the presence of an opiate such as **Heroin**, **Morphine** or **Codeine**.

* This does not exclude the possibility that Fentanyl may be present.



If none of the above colors develop, the above substances are not present in the tested sample.



Disposal

Dilute the acid with some tap water before disposing the kit in the trash.



ENG
v6
01-22



World Leader in
Forensic Products

Manufactured by IDenta Ltd
www.identa-corp.com

Patented Technology • Made in Israel

LOT:

EXP: